

Trocars


Truly customised laparoscopic solutions




A comprehensive range of ergonomically designed trocars, providing simple and safe access into the peritoneal cavity.

- Clear, ribbed cannula for instrument visualisation and effective abdominal wall retention
- Robust universal seal accommodates 5-12mm instruments* whilst maintaining pneumoperitoneum

* Relates to 12mm diameter trocar

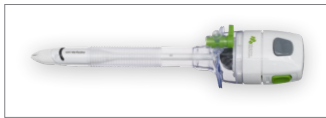

Shielded Bladed Trocar

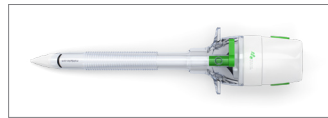


Bladeless Dilating Tip Trocar




Optical Trocar Pistol Grip




Optical Trocar

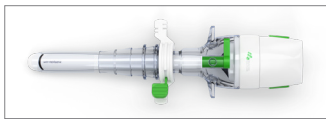


Hasson Trocar

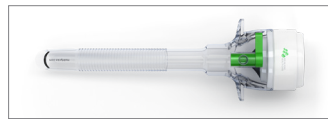

Universal Trocar Cannula

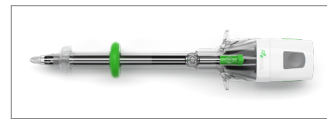

Optical Balloon Trocar

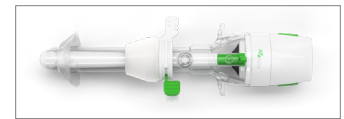

Hasson Balloon Trocar

Shielded bladed trocar

Equipped with a sharp linear blade and a simple intuitive shielding mechanism to prevent inadvertent tissue damage, especially with the 'blind' primary port placement.

- Sharp linear blade allows low insertion force and reduced fascial defect
- Reduced blade width minimises laceration and optimises dilation of tissues

Bladeless trocar

Designed to separate rather than lacerate tissue, causing reduced fascial defect compared with bladed trocars, promoting the reduced requirement for fascial closure.

Optical trocar

Visible trocar allows direct visualisation of the abdominal wall layers as the trocar is inserted, reducing the risk of trocar related injuries especially during primary port placement.

- Provides safe entry into the peritoneal cavity

Hasson trocar

Hasson trocars are designed for surgeons who prefer the security of the open technique for the primary port placement.

- Blunt tip design prevents inadvertent damage to viscera during insertion
- Cannula depth may be rapidly adjusted by means of the simple and effective fixation device

Universal trocar cannulas

Compatible with all same size obturators (shielded bladed, bladeless, optical) providing optimal ordering, storage and cost saving efficiencies.

Available both as sterile single-packed product and for inclusion in Mölnlycke trays.

Balloon Fixation Trocars

- Reduce abdominal intrusion of the cannula which increases surgical working space and optimises visibility, whilst ensuring maintenance of the pneumoperitoneum

Ref. No.	Description	Size	Length	Pcs/Disp.
899300	Shielded Bladed Trocar	5mm	70mm	5
899301	Shielded Bladed Trocar	5mm	100mm	5
899302	Shielded Bladed Trocar	11mm	100mm	5
899303	Shielded Bladed Trocar	11mm	150mm	5
899304	Shielded Bladed Trocar	12mm	100mm	5
899305	Shielded Bladed Trocar	12mm	150mm	5
899308	Bladeless Dilating Tip Trocar	5mm	70mm	5
899309	Bladeless Dilating Tip Trocar	5mm	100mm	5
899310	Bladeless Dilating Tip Trocar	11mm	100mm	5
899311	Bladeless Dilating Tip Trocar	11mm	150mm	5
899312	Bladeless Dilating Tip Trocar	12mm	100mm	5
899313	Bladeless Dilating Tip Trocar	12mm	150mm	5
899336	Bladeless Dilating Tip Trocar	15mm	100mm	5
899314	Optical Trocar Pistol Grip	11mm	100mm	5
899315	Optical Trocar Pistol Grip	12mm	100mm	5
899316	Optical Trocar	5mm	70mm	5
899317	Optical Trocar	5mm	100mm	5

Ref. No.	Description	Size	Length	Pcs/Disp.
899324	Optical Trocar	5mm	150mm	5
899318	Optical Trocar	11mm	100mm	5
899325	Optical Trocar	11mm	150mm	5
899319	Optical Trocar	12mm	100mm	5
899326	Optical Trocar	12mm	150mm	5
899320	Universal Trocar Cannula	5mm	70mm	5
899321	Universal Trocar Cannula	5mm	100mm	5
899322	Universal Trocar Cannula	11mm	100mm	5
899323	Universal Trocar Cannula	12mm	100mm	5
899306	Hasson Trocar	11mm	100mm	5
899307	Hasson Trocar	12mm	100mm	5
899330	Optical Balloon Trocar	5mm	70mm	5
899327	Optical Balloon Trocar	5mm	100mm	5
899331	Optical Balloon Trocar	11mm	100mm	5
899328	Optical Balloon Trocar	12mm	100mm	5
899329	Hasson Balloon Trocar	12mm	100mm	5

Find out more at www.molnlycke.com

Mölnlycke Health Care AB, P.O. Box 13080, Gamlestadsvägen 3 C, SE-402 52 Göteborg, Sweden. Phone + 46 31 722 30 00
The Mölnlycke trademarks, names and logos are registered globally to one or more of the Mölnlycke Health Care Group of Companies. © 2017 Mölnlycke Health Care AB. All rights reserved. HOSU0253

Mölnlycke®